




# BOVIE<sup>®</sup> BANTAM|PRO<sup>™</sup>

## *Electrosurgical Generator*

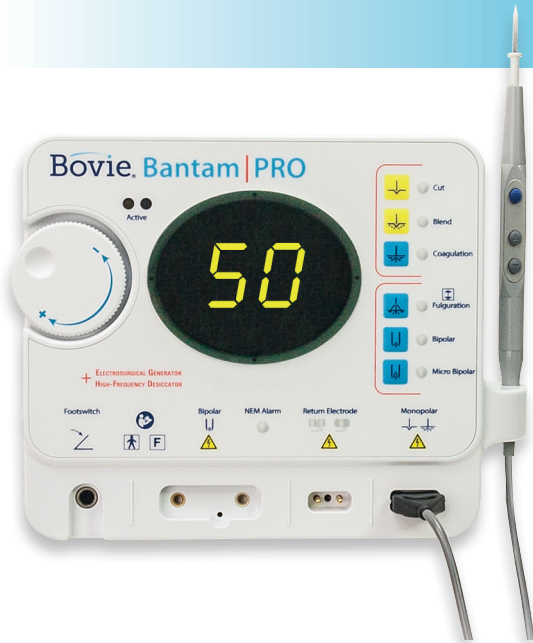

The Bovie Bantam|PRO, A952, is a High Frequency Desiccator with cutting ability. You might call it a split personality. One half of the unit is a high frequency desiccator, an upgrade to our very popular Aaron<sup>®</sup> 940. The other half is a Bovie Generator, giving you 50 watts of cutting power and coag capabilities – plenty of power for most cutting procedure performed in the doctor's office.

High Frequency Desiccator with 50 watts of Cutting capabilities

---

Two machines in one convenient package

---

Adjustable in 1/10th watt increments in the first 10 watts, for the most accurate output available

---

No grounding pad required when working in the desiccation mode

---

When in electrosurgical generator mode, the unit features Bovie NEM<sup>™</sup> pad sensing technology which monitors the return electrode for optimum patient protection.

---



# BOVIE® BANTAM|PRO™ A952 ELECTROSURGICAL GENERATOR



## STANDARD EQUIPMENT

- A952** (1) electrosurgical generator,
- A902** (1) handpiece
- A910** (2) non-sterile disposable handpiece sheaths\*
- A910ST** (2) sterile disposable handpiece sheaths\*
- ES01** (3), **ES02** (1), **ES20** (1) sterile disposable electrodes
- A804** (5), **A806** (5), **A806DE** (2) disposable non-sterile dermal tips\*
- A805** (2), **A807** (2), **A807DE** (2) disposable sterile dermal tips\*
- A1252C** (1) reusable grounding cord
- ESRE** (5) disposable split grounding pads
- A837** wall mount kit\*
- 110 VAC Hospital-grade power cord\*
- User's guide (MC-55-238-001)\*



**A902**



**ES01, ES02, ES20**



**A1252C**



**ESRE**

\* Not Pictured

## COMPLEMENTARY PRODUCTS †



**A813**  
Table Top Stand †



**A803**  
Footswitch †




**A827V**  
Force Cable †



**A824**  
Reusable Bipolar Forceps †



**Lletz Loop**  
Electrodes †

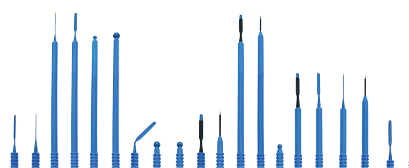

**Resistick II**  
Coated Electrodes †



**Tungsten Wire**  
Loop Electrodes †

† Sold Separately

## PRODUCT SPECIFICATIONS

### DIMENSIONS

Height: 188mm (7.4 in.)  
Width: 228mm (8.98 in.)  
Depth: 105mm (4.13 in.)  
Weight: < 2.26 kg (< 5 lbs.)

### OUTPUT POWER

Cut: 50 watts  
Blend: 50 watts  
Coagulation: 50 watts  
Fulguration: 40 watts  
Bipolar: 40 watts  
Micro Bipolar: 40 watts

### OUTPUT FREQUENCY

304 – 451kHz

### LINE VOLTAGE

A952: 100-240 VAC

### LINE FREQUENCY

50Hz

### WARRANTY

Four Years